## Vulcan® burnout furnaces





Each Vulcan is constructed with a heavy gauge steel cabinet finished with a durable baked enamel coating. The vertical lift door directs heat away from the operator and is standard for safety and space conservation. Each furnace comes with a 1-Year Limited Warranty.



- High performance hybrid muffle
- Longer life, more durable than fiber construction
- Faster heating and cooling than full firebrick
- Integrated door safety switch
- Durable stainless steel front panel
- Digital controller
  - Single-Stage: digital temperature controller (programmable temperature, rate, hold time, alarm, °F/°C operation)
  - Three-Stage: three-stage digital controller (nine multi-stage programs, rate, hold time, delay start, alarm, °F/°C operation)




Available in both single-stage and three-stage options, it has never been easier to select the right burnout furnace for your specific application!

### Technical Data\*

Wide Operating Temperature Range	77-2012°F (25-1100°C)		
Ambient Storage Temperature Range	41-104°F (5+40°C)		
Programmable Heat Rate	1-54°F/min (1-30°C/min)		
Programmable Hold Time	Single-Stage: 0:00-60:00 hrs. / Three-Stage: 0:00-99:59 hrs.		
Muffle Temperature Uniformity	±15°F (±8°C) at steady state		
Operating Temperature (Ambient)	41-104°F (5-40°C)		
Electrical	50-60 Hz		

<sup>\*</sup>These results are based on the testing methods, frequency and procedures of Ransom & Randolph or its approved suppliers. The levels referenced herein are only for general guidance and do not constitute a firm specification.

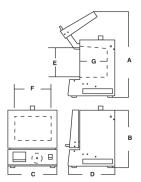

### Exterior Dimensions\*\*

	Model 130	Model 550		
Α	20.0" (510 mm)	25.0" (635 mm)		
В	13.8" (350 mm)	16.0" (410 mm)		
С	11.5" (290 mm)	15.8" (400 mm)		
D	13.0" (330 mm)	17.0" (430 mm)		

<sup>\*\*</sup>Dimension includes exhaust tube.

### Interior Muffle Dimensions

	Model 130	Model 550
Е	4.6" (120 mm)	7.0" (180 mm)
F	5.7" (140 mm)	9.0" (230 mm)
G	5.2" (130 mm)	9.0" (230 mm)



## RANSOM & RANDOLPH

3535 Briarfield Boulevard | Maumee, OH 43537 USA 800.800.7496 | 419.865.9497 | 419.865.9997 (FAX) www.ransom-randolph.com









# Vulcan® burnout furnaces

### Vulcan® Single-Stage Burnout Furnaces

Vulcan single-stage burnout furnaces feature a single-stage controller. The single setpoint digital controller allows the operator to control the rate of temperature climb, and monitor the progress of the program on an easy-to-read, accurate digital display. Adjust the heating rate of the furnace to meet the requirements of a wide range of investments with this easy-to-program option.

Model	<u>S-550</u>		
Item Number	<u>5355730001</u>	5355730002	
Description	S-550 200-240V	S-550 120V	
Volts	200-240 V	100-120 V	
Amps	7.7-9.2 A	16 A	
Watts	1540-2215 W	1920 W	
Watts to Maintain 1832°F (1000°C) (Air Exchange setting at 0)	1050		
Weight	Furnace: 45 lbs. (20 kg) / Shipping: 60 lbs. (27 kg)		

### Vulcan® Three-Stage Burnout Furnaces

Vulcan three-stage burnout furnaces are available in two standard muffle sizes and feature a three-series controller. Threestage digital programmable models have nine multi-stage programs, a hold temperature program, and a delay start option which allows a program's start time to be delayed for up to seven days.

Model	<u>3-130</u>		<u>3-550</u>	
Item Number	<u>5355740001</u>	5355740002	5355740003	5355740004
Description	3-130 200-240V	3-130 100-120V	3-550 200-240V	3-550 100-120V
Volts	200-240V	100-120V	200-240V	100-120V
Amps	4.4 A	12 A	10 A	14.7-16.0 A
Watts	1060 W	1440 W	2400 W	1470-1920 W
Watts to Maintain 1832°F (1000°C) (Air Exchange setting at 0)	525		1050	
Weight	Furnace: 25 lbs. (11 kg) Shipping: 35 lbs. (16 kg)		Furnace: 45 lbs. (20 kg) Shipping: 60 lbs. (27 kg)	

As the conditions or methods of use are beyond our control, we do not assume any responsibility and expressly disclaim any liability for any use of this product. Information contained herein is believed to be true and accurate but all statements or suggestions are made without warranty, expressed or implied, regarding accuracy of the information, the hazards connected with the use of the material or the results to be obtained from the use thereof. Compliance with all applicable federal, state, and local regulations remains the responsibility of the user. All potential liability related to the sale and use of this product is limited to the cost of the particular goods sold in their respective transactions.




RANSOM & RANDOLPH

3535 Briarfield Boulevard | Maumee, OH 43537 USA 800.800.7496 | 419.865.9497 | 419.865.9997 (FAX) www.ransom-randolph.com







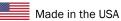

Issue Date: September 22, 2023