R&R Dental



Neycraft™ Multizone™ Duplicating Unit

Quick Start Guide

Before You Begin

Before operating your Neycraft™ MultiZone™ Duplicating Unit, download the R&R Duplicator Connect App. The app will allow you to synchronize your unit with the current date and time and access unit controls.

Note: Ensure Bluetooth is enabled on your mobile device or computer before opening the R&R Duplicator Connect App.

Download the App on Google Play


Download the App on the App Store

Desktop version available under the name R&R Duplicator Connect App in the Microsoft Store.

1 Scan & Select Your Unit


Open the R&R Duplicator Connect App and select **Scan**. When your duplicating unit appears in the app, select it, then select **Info.**

Note: Your unit will appear as RR-1500 followed by a unique seven-character alphanumeric identifier. The identifier will differ from the examples shown below.

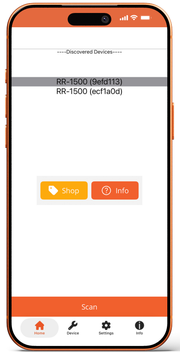


2 Sync Date & Time

From the Info tab, select **Sync RTC.** This synchronizes the unit with the current date and time.

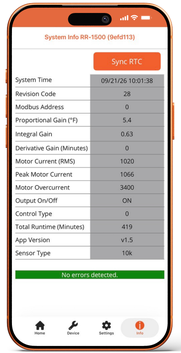


3 Select Your Settings

Choose your preferred unit of measurement. Set the Melt Temp, Idle Temp, and Hold Time according to the duplicating material manufacturer's recommendations, then select **Save.**

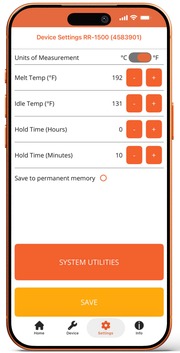

4 Confirm Your Settings

On the **Device** tab, confirm that your settings have saved correctly. Once all settings have been verified, select **Initialize.**

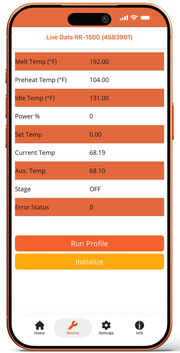

5 Your Unit is Ready

Once the date and time are synchronized and your operating settings have been saved, your unit is ready for operation.



For complete operating instructions, safety information, and R&R duplicating materials compatible with your Neycraft™ MultiZone™ Duplicating Unit, scan the QR code below.

Scan here for duplicating unit operating instructions, safety information, and R&R duplicating materials information.

R&R Ransom&Randolph

Global Headquarters
3535 Briarfield Boulevard
Maumee, Ohio 43537 USA

European Office
Leipziger Straße 40
04571 Rötha Germany

Email
DentalTech@ransom-randolph.com
CustomerService@ransom-randolph.com

Call
US Toll Free: 800-800-7496
US Phone: 419-865-9497
EU: +49 342 06373999

Ransom-Randolph.com